



PRECISION

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one plate to assess.


Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

	Intra-Assay Precision			Inter-Assay Precision		
Sample	1	2	3	1	2	3
n	20	20	20	20	20	20
Mean(ng/ml)	6.113	6.305	6.113	6.314	6.046	6.190
SD	0.052	0.061	0.061	0.069	0.062	0.061
CV(%)	6.899	7.772	7.979	8.817	8.200	7.900

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

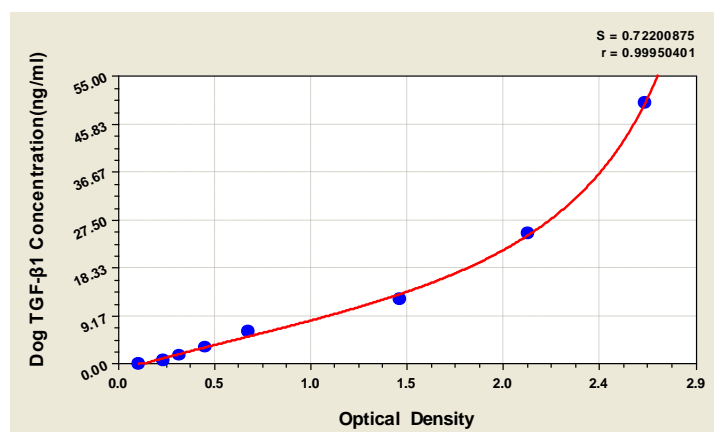

ng/ml	OD1	OD2	Average	corrected
0	0.116	0.112	0.114	
0.781	0.245	0.237	0.241	0.127
1.563	0.331	0.311	0.321	0.207
3.125	0.445	0.458	0.452	0.338
6.25	0.640	0.698	0.669	0.555
12.5	1.436	1.425	1.431	1.317
25	2.123	2.038	2.081	1.967
50	2.701	2.644	2.673	2.559

LOD

0.747 ng/ml

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of dog TGF-β1 in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:1	Average %	97
	Range %	90-109
1:2	Average %	90
	Range %	86-108
1:4	Average %	96
	Range %	91-114
1:8	Average %	92
	Range %	89-99