



PRECISION

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one plate to assess.


Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

Sample	Intra-Assay Precision			Inter-Assay Precision		
	1	2	3	1	2	3
n	20	20	20	20	20	20
Mean(ng/ml)	1.221	1.348	1.289	1.251	1.288	1.288
SD	0.061	0.059	0.065	0.043	0.074	0.067
CV(%)	7.064	6.515	7.331	4.921	8.350	7.560

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

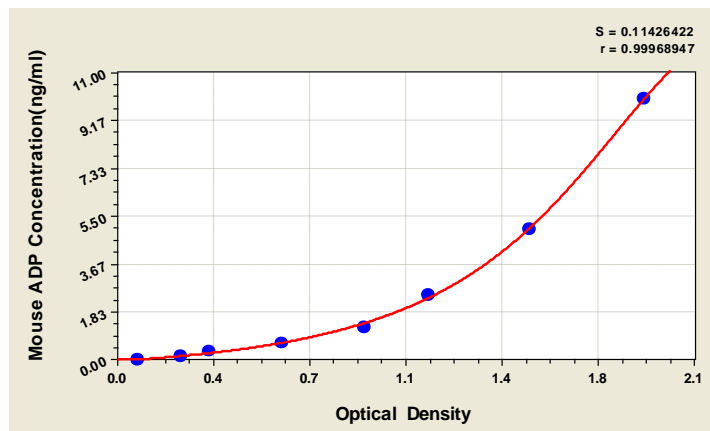

ng/ml	OD1	OD2	Average	corrected
0	0.079	0.089	0.084	
0.156	0.254	0.233	0.244	0.160
0.313	0.341	0.350	0.346	0.262
0.625	0.641	0.592	0.617	0.533
1.25	0.935	0.915	0.925	0.841
2.5	1.167	1.150	1.158	1.074
5	1.540	1.527	1.534	1.450
10	1.990	1.935	1.960	1.876

LOD

0.039 ng/ml

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of mouse ADP in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

Sample	Serum(n=4)
1:4000	Average % 89 Range % 82-104
1:8000	Average % 95 Range % 91-109
1:16000	Average % 92 Range % 85-96
1:32000	Average % 97 Range % 90-112

Sample	Cell culture supernates (n=4)
1:10	Average % 91 Range % 87-103
1:20	Average % 88 Range % 82-109
1:40	Average % 90



1:80

WUHAN HUAMEI BIOTECH CO., LTD.

Postal Address: Wuhan Hi-tech Medical Devices Park, Building B11,
Gaixin Road, Donghu Hi-Tech Development Area, Wuhan, Hubei
Province 430206, P.R. China

Email: cusabio@cusabio.com

Tel: +86-27-87582341

Fax: +86-27-87196150

Range %	84-97
Average %	94
Range %	91-114

