



PRECISION

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one plate to assess.


Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

Sample	Intra-Assay Precision			Inter-Assay Precision		
	1	2	3	1	2	3
n	20	20	20	20	20	20
Mean(ng/ml)	12.155	12.093	11.592	12.324	12.411	12.373
SD	0.033	0.053	0.024	0.034	0.066	0.073
CV(%)	3.946	6.380	2.879	4.104	8.101	8.954

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

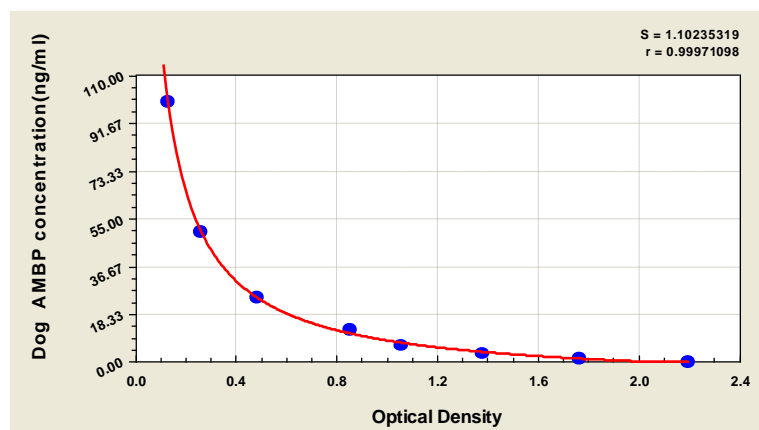

ng/ml	OD1	OD2	Average
100	0.143	0.143	0.143
50	0.288	0.262	0.275
25	0.499	0.505	0.502
12.5	0.901	0.850	0.876
6.25	1.068	1.092	1.080
3.125	1.436	1.376	1.406
1.563	1.815	1.783	1.799
0	2.214	2.254	2.234

LOD

0.649ng/ml

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of Dog AMBP in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:300	Average %	97
	Range %	81-104
1:600	Average %	92
	Range %	81-105
1:1200	Average %	92



1:2400

WUHAN HUAMEI BIOTECH CO., LTD.

Postal Address: Wuhan Hi-tech Medical Devices Park, Building B11,
8818 Gaoxin Road, Donghu Hi-Tech Development Area, Wuhan, Hubei
430206, P.R. China

Email: cusabio@cusabio.com

87-116: +86-27-87582341

Fax: +86-27-87196150

Range %	87-116
Average %	94
Range %	89-102

