


PRECISION
Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one plate to assess.


Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

Sample	Intra-Assay Precision			Inter-Assay Precision		
	1	2	3	1	2	3
n	20	20	20	20	20	20
Mean(ng/ml)	1.234	1.208	1.253	1.225	1.227	1.202
SD	0.026	0.022	0.028	0.046	0.049	0.051
CV(%)	2.377	2.035	2.539	4.224	4.495	4.731

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

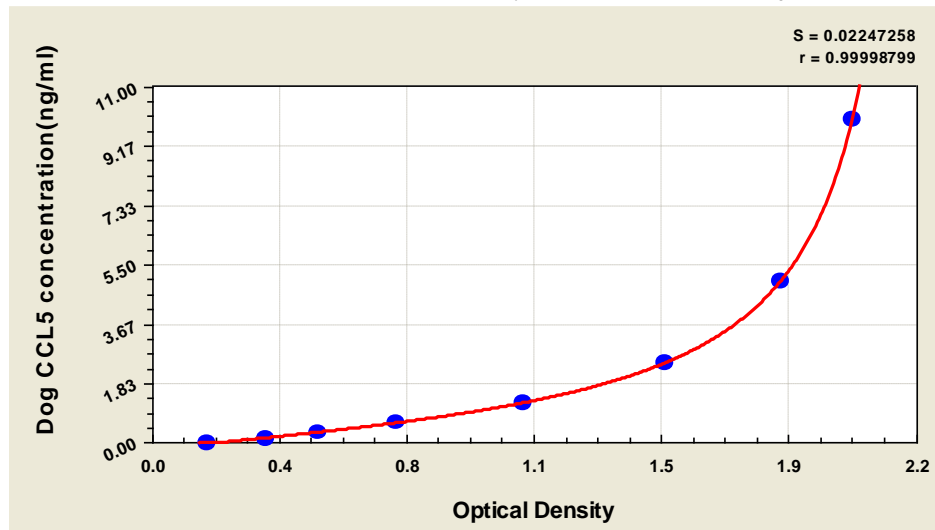

ng/ml	OD1	OD2	Average	corrected
0	0.184	0.165	0.175	
0.156	0.346	0.347	0.347	0.173
0.313	0.520	0.473	0.497	0.322
0.625	0.715	0.729	0.722	0.548
1.25	1.075	1.105	1.090	0.916
2.5	1.507	1.493	1.500	1.326
5	1.899	1.771	1.835	1.661
10	2.049	2.032	2.041	1.866

LOD

0.081 ng/ml

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of Dog CCL5 various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:2	Average %	83
	Range %	80-94
1:4	Average %	85
	Range %	82-99
1:8	Average %	92
	Range %	84-101
1:16	Average %	95
	Range %	89-105